

sCUBE[®] RXN, P&F Module




F-18 Reaction,
Purification & Formulation



Registered ISO 9001:2008 EN/ISO 13485:2012

Compact, Built-in HPLC System

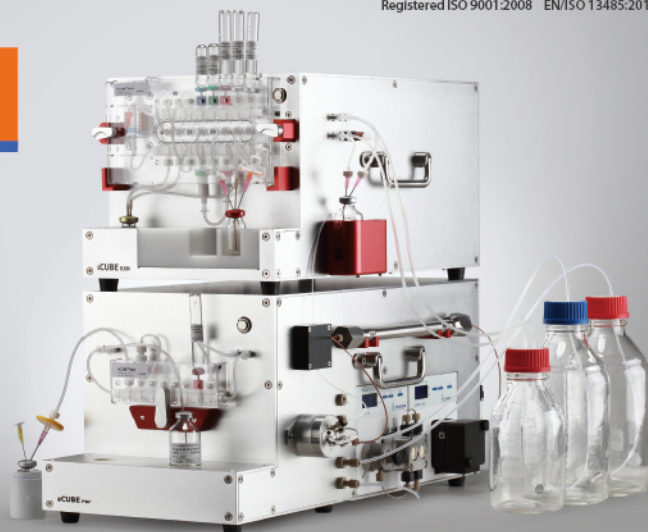

Commercial
Production
reliability



Compact
Size



EASY
"Disposable
Cassette type"



GMP
Compliance



Multi-tracer
application


sCUBE RXN, P&F system is designed to

“^{[18]F}F- to ^{[18]F}Tracer Tracking”


8 - Radiation Detectors(RI)

- sCUBE RXN module(4; QMA, SPE, Reactor, Vial)
- sCUBE P&F module(4; SPE, INJECT, Final Vial, HPLC)

“Built-in HPLC”, auto-Injection system

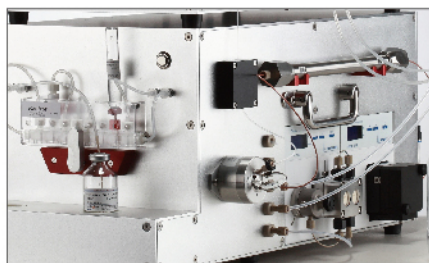

Load complete

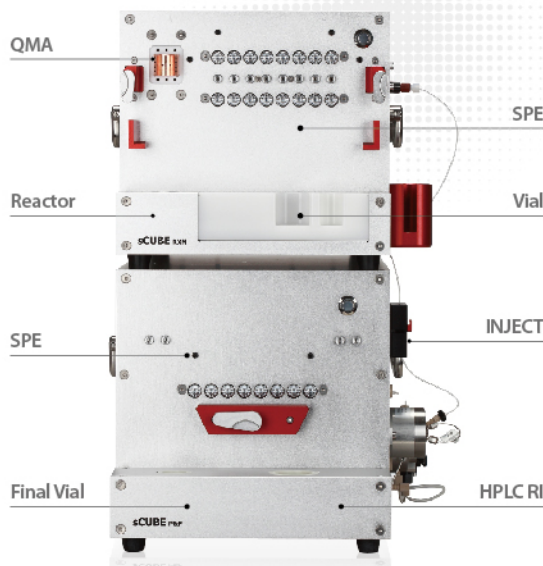


sCUBE® RXN, P&F Module

Specifications




1. Disposable cassette, Reagent	- Multi-tracer possible - Easily to set up and remove
2. GMP compliance	Suitability of GMP facility (IQ, OQ, PQ)
3. Compact size	250(W) x 476(H) x 550(D)mm
4. Built-in HPLC system	High specific activity
5. Auto Injection system	Electrical Auto Injector with RI Detector
6. Commercial production reliability	[¹⁸F]FP-CIT (for Parkinson's Disease imaging) Perform Commercial production for 3 years 24.3% ± 2.5 (50 min, n.d.c.), n=821, ~ Ci scale/day [¹⁸F]FC119S (New Tracer) Perform Routine production for 2 years (Clinical trial) 27.7% ± 4 (55 min, n.d.c.), n=70, ~ Ci scale/day
7. Multi-tracer application	[¹⁸F]FLT (Commercial production) 22.2% ± 6.4 (55 min, n.d.c.), n=24, ~ Ci scale/day [¹⁸F]FDG (using sCUBE RXN module only) ~ 60% (30 min, n.d.c.)

Performance of sCUBE RXN, P&F module

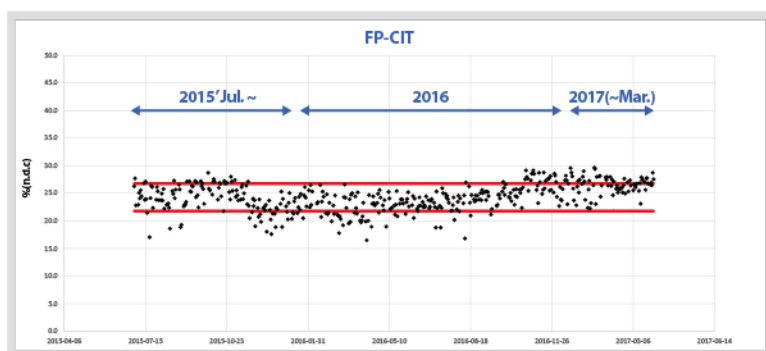

[¹⁸F]FP-CIT(PDvue®)
N=821, 24.3±2.5%
Total Product(FP-CIT) : 413 Ci

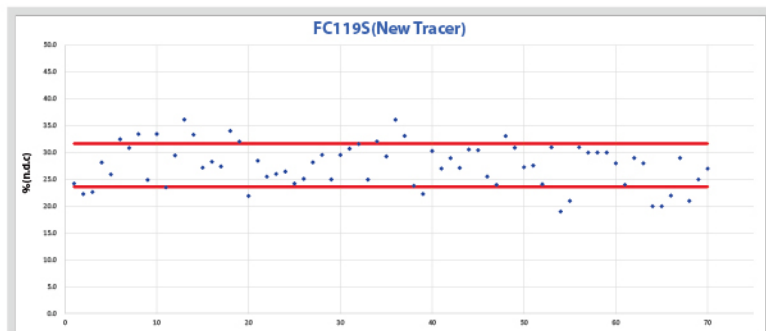

[¹⁸F]FC119S(new tracer)
N=70, 27.7±4%
Total Product(FC119S) : 9.2 Ci